



# SAFETY STORAGE CABINETS



# CONTENTS

<b>Safety Storage Cabinet Features</b>	<b>3</b>
--	----------

## **Model**

Flammable Safety Cabinets	5
Corrosive Safety Cabinets	6
Combustible Safety Cabinets	7
Hazardous Safety Cabinets	8
Polypropylene Safety Cabinets	9
Stainless Steel Safety Cabinets	10
Flammable Drum Storage Cabinets	11

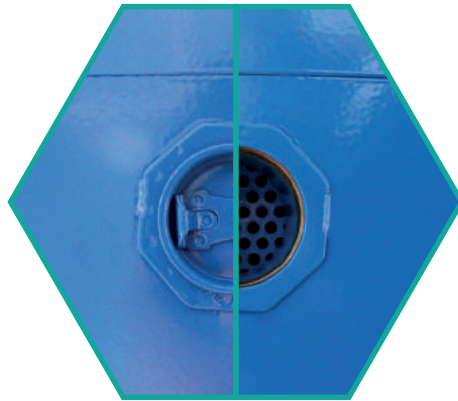
<b>Warranty</b>	<b>12</b>
-----------------	-----------

# Safety Storage Cabinet Features



- Constructed from industrial grade electro-galvanised steel 18 gauge for strength and durability
- Epoxy powder coated to meet harsh industrial/laboratory environments
- Double layered construction with 1.5" airspace between layers
- Fully welded 2" sill to contain chemical spills within the cabinet
- 2 ventilation ports protected with flash arrestors
- 3 point door latching mechanism
- Large Warning Labels for easy identification
- Adjustable shelves with 1/2" centers
- Complies and 3rd party lab tested to NFPA 30 6.3.3 para 2

**Adjustable galvanised steel shelves with heavy duty stainless steel pins**




**2" vents with flash arrestors (exterior)**


**Double wall construction**



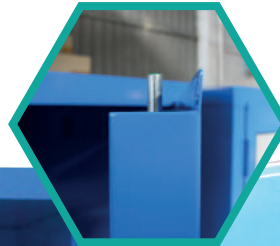
**3 pin lock**



**2" deep spill sump**

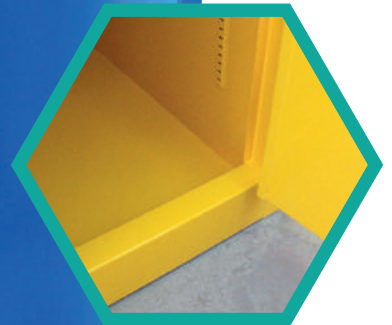

**Heavy duty hinges**

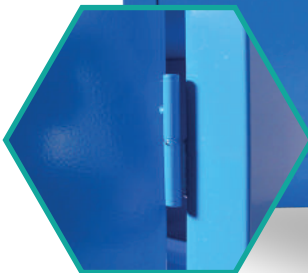


**3 point door latch**









**Optional PP tray**

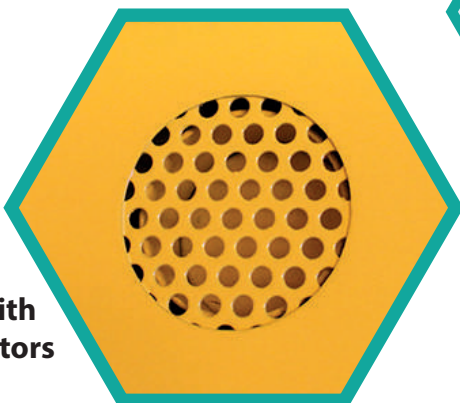**2" vents with flash arrestors (interior)**








# Flammable Safety Cabinets

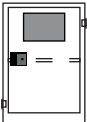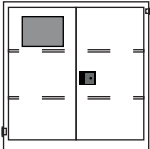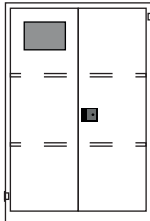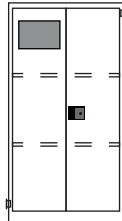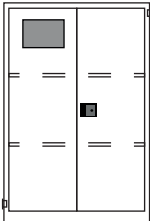
MODEL	AD15F	AD30F	AD45F	AD60F	AD90F
Gallon	15	30	45	60	90
Dimensions, mm (LxWxH)	 610mm x 550mm x 908mm	 1092mm x 457mm x 1118mm	 1092mm x 457mm x 1651mm	 864mm x 864mm x 1651mm	 1092mm x 864mm x 1651mm
Weight	75 Kg	105 Kg	145 Kg	180 Kg	220 Kg
Design	Manual closing door				
Door	1 door	2 door	2 door	2 door	2 door
Shelves	1	2	2	2	2

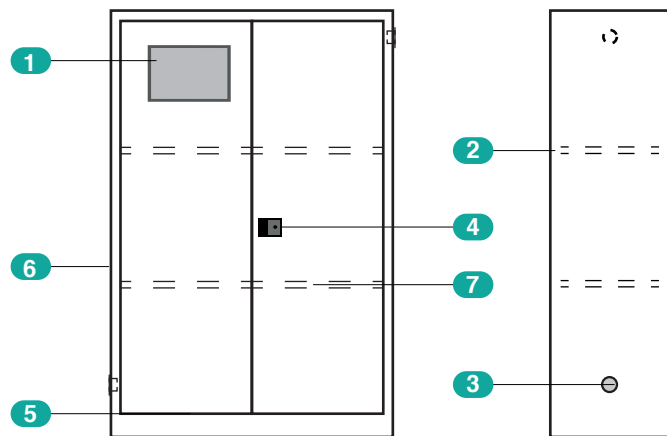

1. Safety Warning Sticker
2. Double wall construction with 1-1/2" air space
3. 2 nos 2" vents with flash arrestors
4. Flush mounted locking handle
5. Leak proof sill 2" deep, prevents leakage in the event of accidental spills
6. Epoxy powder coated
7. Galvanised steel shelves adjusted on 1/2" centers






# Corrosive Safety Cabinets

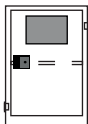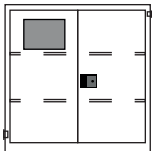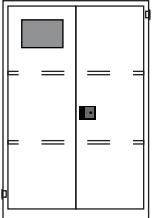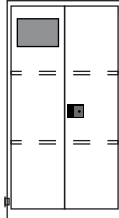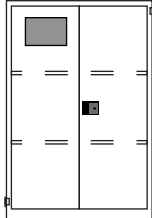
MODEL	AD15C	AD30C	AD45C	AD60C	AD90C
Gallon	15	30	45	60	90
Dimensions, mm (LxWxH)	 610mm x 550mm x 908mm	 1092mm x 457mm x 1118mm	 1092mm x 457mm x 1651mm	 864mm x 864mm x 1651mm	 1092mm x 864mm x 1651mm
Weight	75 Kg	105 Kg	145 Kg	180 Kg	220 Kg
Design	Manual closing door				
Door	1 door	2 door	2 door	2 door	2 door
Shelves	1	2	2	2	2

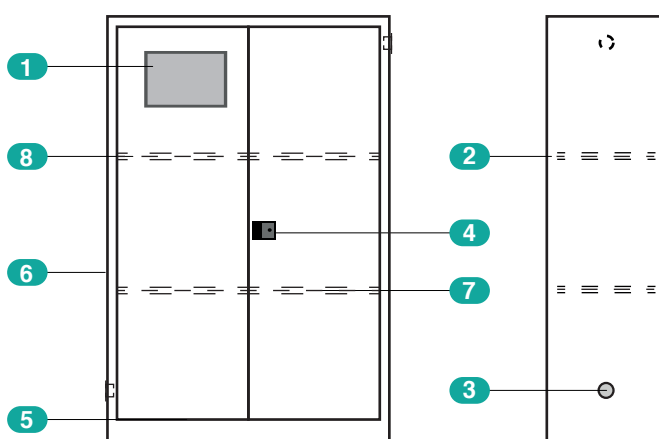

1. Safety Warning Sticker
2. Double wall construction with 1-1/2" air space
3. 2 nos 2" vents with flash arrestors
4. Flush mounted locking handle
5. Leak proof sill 2" deep, prevents leakage in the event of accidental spills
6. Epoxy powder coated
7. Galvanized steel shelves adjusted on 1/2" centers
8. PP Trays






# Combustible Safety Cabinets

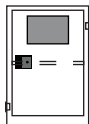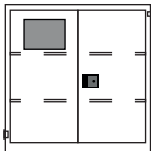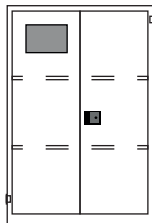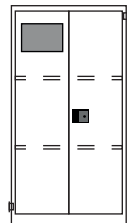
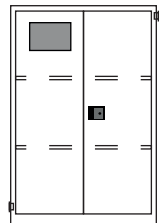
MODEL	AD15R	AD30R	AD45R	AD60R	AD90R
Gallon	15	30	45	60	90
Dimensions, mm (LxWxH)	 610mm x 550mm x 908mm	 1092mm x 457mm x 1118mm	 1092mm x 457mm x 1651mm	 864mm x 864mm x 1651mm	 1092mm x 864mm x 1651mm
Weight	75 Kg	105 Kg	145 Kg	180 Kg	220 Kg
Design	Manual closing door				
Door	1 door	2 door	2 door	2 door	2 door
Shelves	1	2	2	2	2

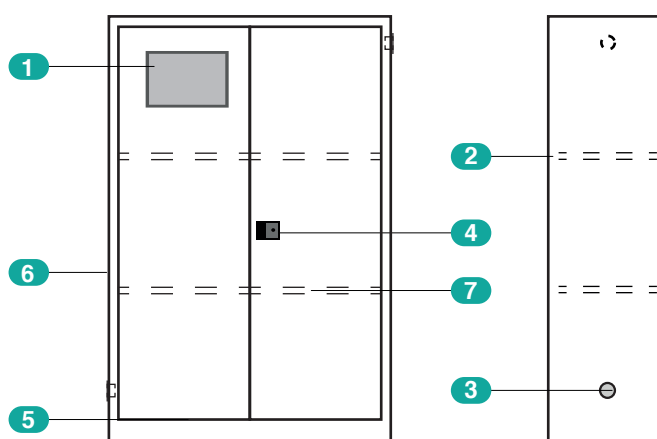

1. Safety Warning Sticker
2. Double wall construction with 1-1/2" air space
3. 2 nos 2" vents with flash arrestors
4. Flush mounted locking handle
5. Leak proof sill 2" deep, prevents leakage in the event of accidental spills
6. Epoxy powder coated
7. Galvanised steel shelves adjusted on 1/2" centers

# Hazardous Safety Cabinets

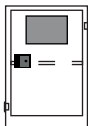
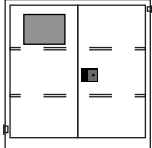
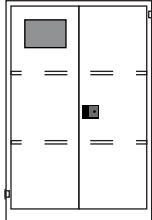
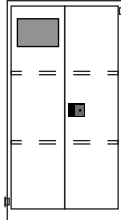
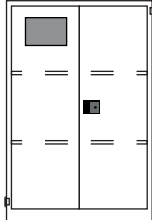
MODEL	AD15H	AD30H	AD45H	AD60H	AD90H
Gallon	15	30	45	60	90
Dimensions, mm (LxWxH)	 610mm x 550mm x 908mm	 1092mm x 457mm x 1118mm	 1092mm x 457mm x 1651mm	 864mm x 864mm x 1651mm	 1092mm x 864mm x 1651mm
Weight	75 Kg	105 Kg	145 Kg	180 Kg	220 Kg
Design	Manual closing door				
Door	1 door	2 door	2 door	2 door	2 door
Shelves	1	2	2	2	2

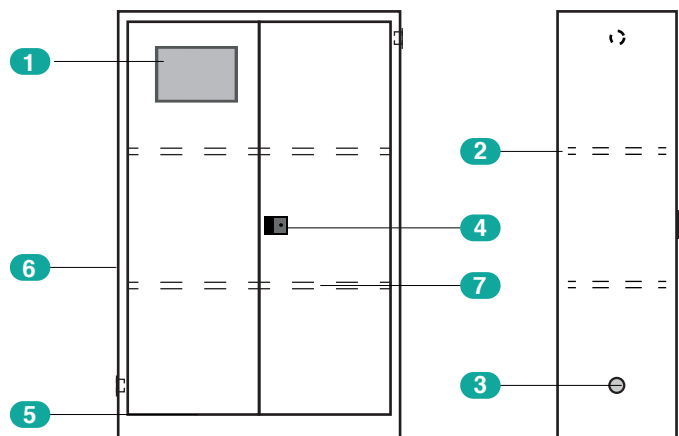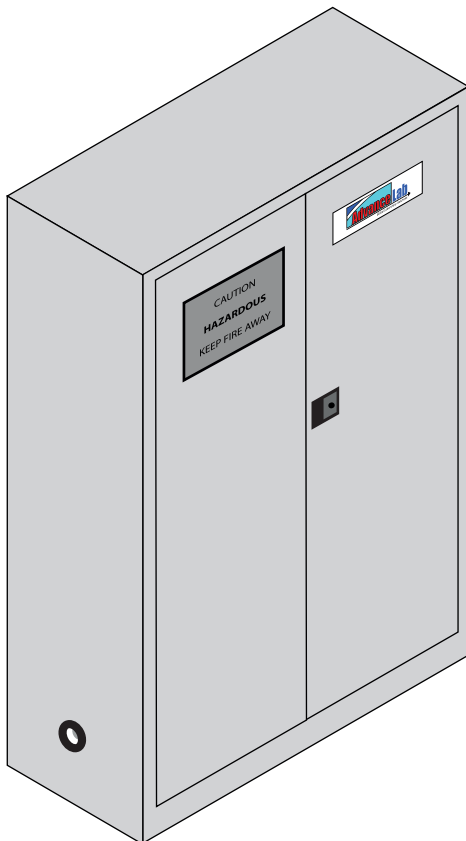


1. Safety Warning Sticker
2. Double wall construction with 1-1/2" air space
3. 2 nos 2" vents with flash arrestors
4. Flush mounted locking handle
5. Leak proof sill 2" deep, prevents leakage in the event of accidental spills
6. Epoxy powder coated
7. Galvanised steel shelves adjusted on 1/2" centers



# Polypropylene Safety Cabinets

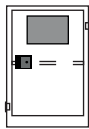
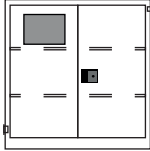
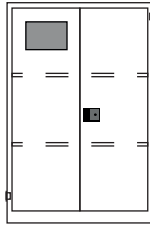
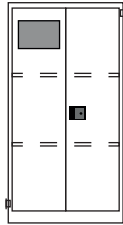
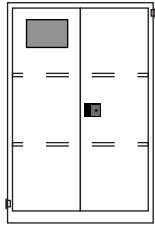
MODEL	AD15PP	AD30PP	AD45PP	AD60PP	AD90PP
Gallon	15	30	45	60	90
Dimensions, mm (LxWxH)	 610mm x 550mm x 908mm	 1092mm x 457mm x 1118mm	 1092mm x 457mm x 1651mm	 864mm x 864mm x 1651mm	 1092mm x 864mm x 1651mm
Weight	75 Kg	105 Kg	145 Kg	180 Kg	220 Kg
Design	Manual closing door				
Door	1 door	2 door	2 door	2 door	2 door
Shelves	1	2	2	2	2

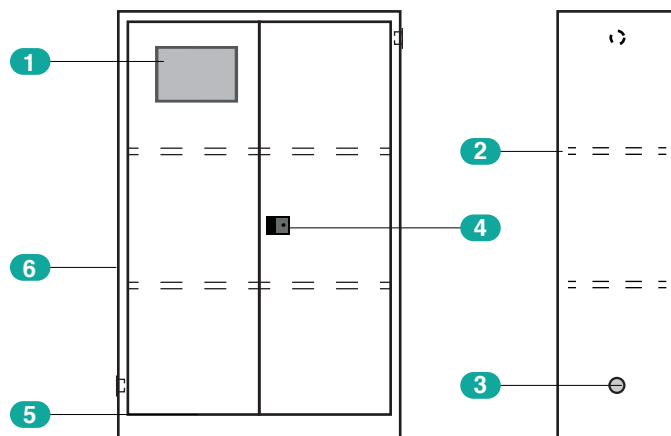


1. Safety Warning Sticker
2. Full PP hot weld construction (Single Layer)
3. 2 nos 2" vents with flash arrestors
4. Flush mounted locking handle
5. Leak proof sill 2" deep, prevents leakage in the event of accidental spills
6. PP adjustable shelves

\* Single layer and not compliant to NFPA 30




# Stainless Steel Safety Cabinets

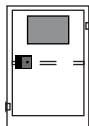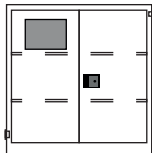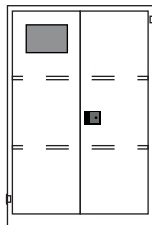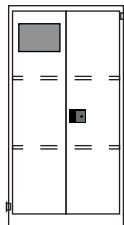
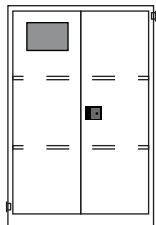
MODEL	AD15SS	AD30SS	AD45SS	AD60SS	AD90SS
Gallon	15	30	45	60	90
Dimensions, mm (LxWxH)	 610mm x 550mm x 908mm	 1092mm x 457mm x 1118mm	 1092mm x 457mm x 1651mm	 864mm x 864mm x 1651mm	 1092mm x 864mm x 1651mm
Weight	75 Kg	105 Kg	145 Kg	180 Kg	220 Kg
Design	Manual closing door				
Door	1 door	2 door	2 door	2 door	2 door
Shelves	1	2	2	2	2

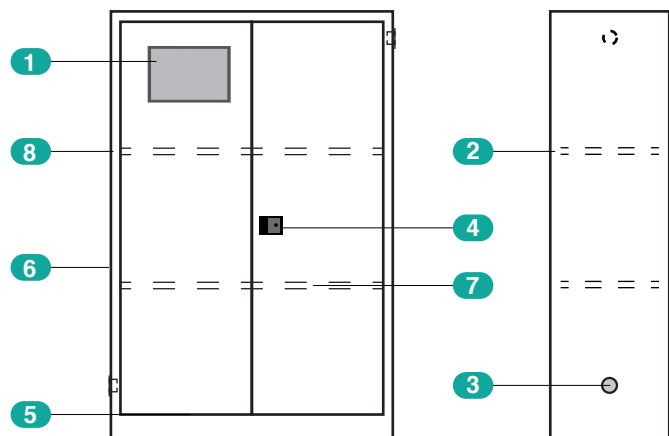

1. Safety Warning Sticker
2. Double wall construction with 1-1/2" air space
3. 2 nos 2" vents with flash arrestors
4. Flush mounted locking handle
5. Leak proof sill 2" deep, prevents leakage in the event of accidental spills
6. Stainless steel construction
7. Stainless steel shelves adjusted on 1/2" centers
8. Constructed from SS304 with 4B finish (SS316 Optional)

# Flammable Drum Storage Cabinets

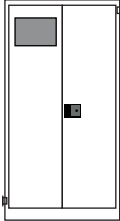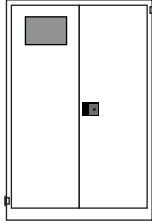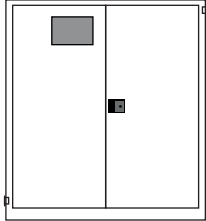MODEL	AD55DS1	AD30DS2	AD55DS2
Dimensions, mm (LxWxH)	864 x 864 x 1651 	1092 x 864 x 1651 	1496 x 864 x 1651 
Weight	150 Kg	220 Kg	280 Kg
Design	Manual closing door		
Door	2 door	2 door	2 door
Capacity	1 x 55 gallon drum	2 x 30 gallon drum	2 x 55 gallon drum

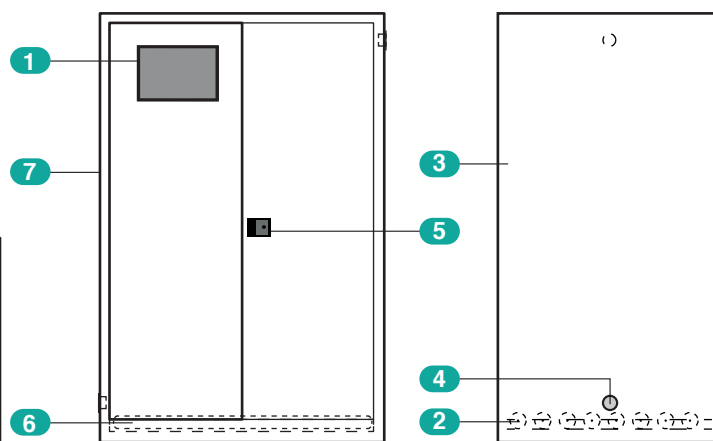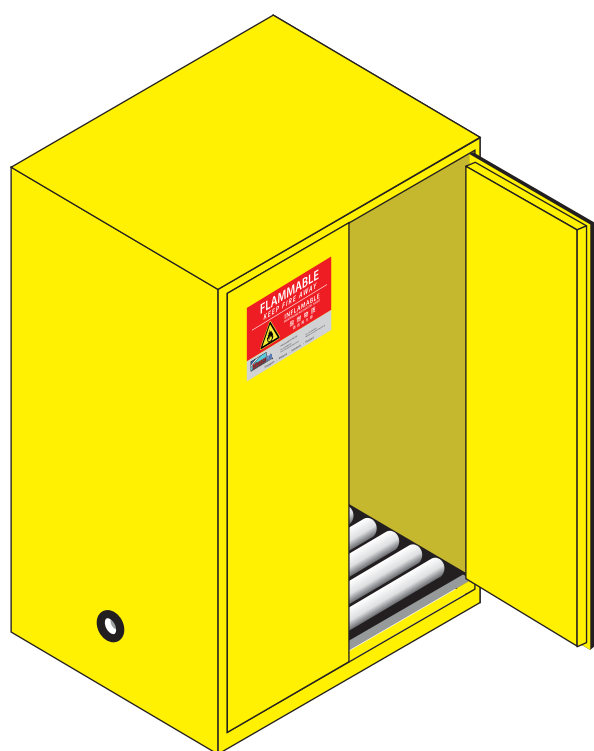

1. Safety Warning Sticker
2. 60mm dia rollers, steel
3. Double wall construction with 1-1/2" air space
4. 2 nos 2" vents with flash arrestors
5. Flush mounted locking handle
6. Leak proof sill 2" deep, prevents leakage in the event of accidental spills
7. Epoxy powder coated

# Warranty

**Products Manufactured by Advancelab:** Advancelab(S) Pte.Ltd., warrants products that it manufactures to be free from defects for a period of 12 months for parts, commencing from the date of shipment. Advancelab's sole responsibility is to repair or replace, at its option, any part of the product that proves defective or malfunctioning during this time limit. This warranty is void if the equipment is abused or modified by the customer, is operated outside Advancelab's operating instructions or specifications, or is used in any application other than that for which it is specified. This warranty does not include routine maintenance or service procedures, breakage, shipping damage, nor damage from misuse, intentional or unintentional abuse, neglect, natural disasters, or acts of God.

**Freight Shortage or Damage:** Upon receipt of any equipment from Advancelab, customer shall immediately unpack and inspect for damage or shortage. The customer shall not accept a damaged package or a short shipment until the carrier makes a "damage or shortage" notation on both the carrier's and customer's copy of the freight bill or delivery receipt. Service title passes when the shipment is loaded, so customer is responsible for filing and collecting a freight claim. Any replacement products must be ordered and paid for separately.

Generally, customers can improve the chance of collecting on a freight claim by following these procedures:

1. Formally requesting that the carrier inspect the shipment immediately upon suspecting damage or shortage to verify condition.
2. Notifying the carrier upon discovery of concealed damage and requesting an inspection within 15 days of receipt, both in person or phone and following up via mail.
3. Keeping the shipment as intact as possible, including retaining original packaging materials and keeping the product as close to the original receiving location as possible.
4. Holding salvage for disposition by the carrier.

**All Claims:** Advancelab (S) Pte.Ltd., expressly disclaims all other warranties, expressed or implied or implied by statute, including the warranties of merchantability or fitness for intended use. Advancelab is not responsible for consequential or incidental damages arising out of the purchase or use of the products supplied by Advancelab. Advancelab is not liable for damage to facilities, other equipment, products, property or personnel of others, or of their agents, suppliers, or affiliated parties, which is caused or alleged to have been caused by products supplied by Advancelab. In any event or series of events, Advancelab's total liability for any and all damages whatsoever is limited to the lesser of the actual damages or the original invoice cost of the items alleged to have caused the damage. The customer's sole and exclusive remedy for any cause of action whatsoever is repair or replacement of the non-conforming products or refund of the actual purchase price, at the sole option of Advancelab. All claims must be made in writing within 90 days of the date the product was shipped. Any claims not made within this time limit shall be deemed waived by the customer. Advancelab is not responsible for any additional costs of repair caused by poor packaging or in-shipment damage during return.

**Warranty Returns:** All warranty returns must be authorized in advance by Advancelab and approved by writing. Unless approved in advance for good reason, all returns must be in original condition, including all manuals, and must be packaged in original packaging materials. All returned goods are to be shipped to Advancelab, freight prepaid at customer's expense.



# SAFETY STORAGE RACKS CABINETS

## Singapore

Advancelab (S) Pte Ltd

No.52 Senang Crescent,  
Singapore 416619  
Tel: +65 6448 8255  
Fax: +65 6448 9833  
Email: [info@advancelab.com.sg](mailto:info@advancelab.com.sg)  
Website: [www.advancelab.com.sg](http://www.advancelab.com.sg)

## Malaysia

Advancelab Sdn Bhd

KL Office: Unit B-1-10, Block B, No. 2,  
Jalan PJU 1A/7A, Oasis Square, Ara  
Damansara, 47301 Petaling Jaya,  
Selangor Darul Ehsan, Malaysia.  
Tel: +603 7831 0188  
Fax: +603 78310588  
Email: [info-my@advancelab-global.com](mailto:info-my@advancelab-global.com)  
Website: [www.advancelab.com.sg](http://www.advancelab.com.sg)

Johor Factory: No. 3388, Jalan  
Pekeliling Tanjung 27/2, Kawasan  
Perindustrian Indahpura, 81000  
Kulaijaya, Johor, Malaysia.  
Tel: +607 660 8877  
Fax: +607 660 8866  
Email: [info-my@advancelab-global.com](mailto:info-my@advancelab-global.com)  
Website: [www.advancelab.com.sg](http://www.advancelab.com.sg)

## Myanmar

Advancelab Scientific &  
Engineering Co., Ltd

No(81/2), 7th Street, Than Thu Mar Road,  
(14/1)Quarter, South Okkalapa Township,  
11091 Yangon, Myanmar.  
Tel: +95 (1) 572393 / +95 (9) 779753802  
Fax: +95 (1) 572393  
Email: [info-mm@advancelab-global.com](mailto:info-mm@advancelab-global.com)  
Website: [www.advancelab.com.sg](http://www.advancelab.com.sg)

## Indonesia

PT. Advancelab Saintifik

Jl. Boulevard Raya Blok QJ 5 No. 23  
Kelapa Gading, Jakarta Utara, 14240  
Indonesia  
Tel: +6221 45854569 / +6221 45854570  
Fax: +6221 45854570  
Email: [info-id@advancelab-global.com](mailto:info-id@advancelab-global.com)  
Website: [www.advancelab.com.sg](http://www.advancelab.com.sg)

## Thailand

Advancelab (Thailand)  
Co., Ltd.

147 / 170 - 171 Baromrajchonnanee  
Road, Arunamarin Bangkoknoi, Bangkok  
10700, Thailand.  
Tel: +66 2 884 9480  
Fax: +66 2 884 6441  
Email: [primasci@advancelab-global.com](mailto:primasci@advancelab-global.com)  
Website: [www.advancelab.com.sg](http://www.advancelab.com.sg)

## U.A.E (Dubai)

Advancelab FZCO

No.FZJOAB1310  
Jafza One, Jebel Ali Dubai  
United Arab Emirates  
Email: [info@advancelab.com.sg](mailto:info@advancelab.com.sg)  
Website: [www.advancelab.com.sg](http://www.advancelab.com.sg)

## Vietnam

Advancelab Pte Ltd

83/10 Bach Dang Street, Ward 2,  
Tan Binh District,  
Ho Chi Minh City 700000, Vietnam  
Tel: +84 8 2244 6300  
Fax: +84 8 3845 1974  
Email: [info@advancelab.com.sg](mailto:info@advancelab.com.sg)  
Website: [www.advancelab.com.sg](http://www.advancelab.com.sg)

